



Standard Test Method for Quantification of a *Pseudomonas aeruginosa* Biofilm Grown with Shear and Continuous Flow Using a Rotating Disk Reactor¹

This standard is issued under the fixed designation E 2196; the number immediately following the designation indicates the year of original adoption or, in the case of revision, the year of last revision. A number in parentheses indicates the year of last reapproval. A superscript epsilon (ε) indicates an editorial change since the last revision or reapproval.

1. Scope

1.1 This test method is used for growing a repeatable² *Pseudomonas aeruginosa* biofilm in a continuously stirred flow reactor. In addition, the test method describes how to sample and analyze biofilm for viable cells.

1.2 In this test method, biofilm population density is recorded as log colony forming units per surface area.

1.3 Basic microbiology training is required to perform this test method. *This standard does not claim to address all of the safety problems associated with its use. It is the responsibility of the user of this standard to establish appropriate safety practices and determine the applicability of regulatory limitations prior to use.*

2. Referenced Documents

2.1 Other Standards:

Buffered Dilution Water Preparation —Method 9050 C.1a³
Rotating Disk Reactor —Repeatability and Relevance⁴
Rotating Disk Reactor —Efficacy Test Method⁵

3. Terminology

3.1 *biofilm, n*—microorganisms living in a self-organized, cooperative community attached to surfaces, interfaces, or each other, embedded in a matrix of extracellular polymeric sub-

stances of microbial origin, while exhibiting an altered phenotype with respect to growth rate and gene transcription.

3.1.1 *Discussion*—Biofilms may be comprised of bacteria, fungi, algae, protozoa, viruses, or infinite combinations of these microorganisms. The qualitative characteristics of a biofilm, including, but not limited to, population density, taxonomic diversity, thickness, chemical gradients, chemical composition, consistency, and other materials in the matrix that are not produced by the biofilm microorganisms; are controlled by the physiochemical environment in which it exists.

3.2 *coupon*—biofilm sample surface.

4. Summary of Test Method

4.1 This test method is used for growing a repeatable *Pseudomonas aeruginosa* biofilm in a rotating disk reactor. The biofilm is established by operating the reactor in batch mode (no flow) for 24 h. A steady state growth (attachment is equal to detachment) is reached while the reactor operates for an additional 24 h with continuous flow of the nutrients. The residence time of the nutrients in the reactor is set to select for biofilm growth, and is species and reactor parameter specific. During the entire 48 h, the biofilm is exposed to continuous fluid shear from the rotation of the disk. At the end of the 48 h, biofilm accumulation is quantified by removing coupons from the disk, scraping the biofilm from the coupon surface, disaggregating the clumps, then diluting and plating for viable cell enumeration.

5. Significance and Use

5.1 Bacteria that exist in a biofilm are phenotypically different from suspended cells of the same genotype. The study of biofilm in the laboratory requires protocols that account for this difference. Laboratory biofilms are engineered in growth reactors designed to produce a specific biofilm type. Altering system parameters will correspondingly result in a change in the biofilm. The purpose of this method is to direct a user in the laboratory study of biofilms by clearly defining each system parameter. This method will enable a person to grow, sample, and analyze a laboratory biofilm.

¹ This test method is under the jurisdiction of ASTM Committee E35 on Pesticides and Alternative Control Agents and is the direct responsibility of Subcommittee E35.15 on Antimicrobial Agents.

Current edition approved April 1, 2007. Published May 2007. Originally approved in 2002. Last previous edition approved in 2002 as E 2196 – 02.

² Ellison, S.L.R., M. Rosslein, A. Williams. (Eds.) 2000. *Quantifying Uncertainty in Analytical Measurement*, 2nd Edition. Eurachem.

³ Eaton, A.D., L.S. Clesceri, Rice, E.W., A.E. Greenberg. (Eds.) *Standard Methods for the Examination of Water and Waste Water*, 21st Edition. American Public Health Association, American Water Works Association, Water Environment Federation. Washington D.C., 2005.

⁴ Zelter, N., M. Hamilton, B. Pitts, D. Goeres, D. Walker, P. Sturman, J. Heersink. 1999. Methods for measuring antimicrobial effects on biofilm bacteria: from laboratory to field. In: Doyle, R.J. (Ed.), *Methods in Enzymology-Biofilms* Vol 310, Academic Press, San Diego, CA, pp. 608-628.

⁵ Zelter, N., M. Hamilton, D. Goeres, J. Heersink. 2001. Development of a Standardized Antibiofilm Test. In: Doyle, R.J. (Ed.), *Methods in Enzymology-Biofilms* Vol 337, Academic Press, San Diego, CA, pp. 363-376.

6. Apparatus

- 6.1 *Wooden Applicator Sticks*, sterile.
 - 6.2 *Inoculating Loop*.
 - 6.3 *Petri Dish*, 100 by 15 mm, plastic, sterile and empty to hold rotor while sampling.
 - 6.4 *Culture Tubes and Culture Tube Closures*, any with a volume capability of 10 mL and diameter no less than 6 cm. Recommended size is 16 by 125 mm borosilicate glass with threaded opening.
 - 6.5 *Pipetter*, continuously adjustable pipetter with volume capability of 1 mL.
 - 6.6 *Vortex*, any vortex that will ensure proper agitation and mixing of culture tubes.
 - 6.7 *Homogenizer*, any capable of mixing at $20\,500 \pm 5000$ r/min in a 5 to 10 mL volume.
 - 6.8 *Homogenizer Probe*, any capable of mixing at $20\,500 \pm 5000$ r/min in a 5 to 10 mL volume and can withstand autoclaving or other means of sterilization.
 - 6.9 *Sonicator*, any noncavitating sonicating bath that operates at 50 to 60 Hz.
 - 6.10 *Syringe*, sterile, 1 mL syringe used during reactor inoculation.
 - 6.10.1 *Needle*, sterile, 22 gauge needle used with syringe to inoculate reactor.
 - 6.11 *Bunsen Burner*, used to flame inoculating loop and other instruments.
 - 6.12 *Stainless Steel Dissecting Tools*.
 - 6.13 *Stainless Steel Hemostat Clamp with Curved Tip*.
 - 6.14 *Environmental Shaker*, capable of maintaining temperature of $35 \pm 2^\circ\text{C}$.
 - 6.15 *Analytical Balance*, sensitive to 0.01 g.
 - 6.16 *Sterilizers*, any steam sterilizer capable of producing the conditions of sterilization is acceptable.
 - 6.17 *Colony Counter*, any one of several types may be used, such as the Quebec, Buck, and Wolfhugel. A hand tally for the recording of the bacterial count is recommended if manual counting is done.
 - 6.18 *Peristaltic Pump*, pump head capable of holding tubing with ID 3.1 mm and OD 3.2 mm.
 - 6.19 *Magnetic Stir Plate*, top plate 10.16 by 10.16 cm, capable of rotating at 200 ± 100 r/min.
- NOTE 1—R/min may be measured using a strobe light.
- 6.20 *Silicone Tubing*, two sizes of tubing: one with ID 3.1 mm and OD 3.2 mm and the other with ID 7.9 mm and OD 9.5 mm. Both sizes must withstand sterilization.
 - 6.21 *Glass Flow Break*, any that will connect with tubing of ID 3.1 mm and withstand sterilization.
 - 6.21.1 *Clamp*, used to hold flow break, extension clamp with 0.5 cm minimum grip size.
 - 6.21.2 *Clamp Stand*, height no less than 76.2 cm, used with clamp to suspend glass flow break vertically and stabilize tubing above reactor.
 - 6.22 *Reactor Components*⁶:

6.22.1 *Berzelius Pyrex Beaker*, 1000 mL without pour spout, 9.5 ± 0.5 cm diameter. Pyrex barbed outlet spout added at 250 mL $\pm$ 15 mL mark at 30 to 45° angle, spout should accommodate silicone tubing with an ID of 8 to 11 mm.

NOTE 2—The rotor, described in 6.22.3, will displace approximately 50 mL of liquid. Therefore, an outlet spout at the 250 mL mark will result in approximately a 200 mL operating volume. The user is encouraged to confirm the actual liquid volume in the reactor, when the rotor is in place, before use. The measured volume is used to calculate an exact pump flow rate.

6.22.2 *Reactor Top*, size 15 rubber or machined stopper, 3 to 4 holes bored through stopper to accommodate 6 cm pieces of fire-polished glass tubing or other suitable rigid autoclavable tubing with OD 4 to 6 mm, as shown in Fig. 1. Another hole can be added to the stopper to contain an inoculum port. The inoculum port consists of a 6 cm piece of fire-polished glass tubing or other suitable rigid autoclavable tubing fitted with a septum.

6.22.3 *Rotor or Disk*, nominal (see Note 3) 1.6 mm thick PTFE sheet cut into a disk with a diameter of 7.0 ± 0.2 cm containing 6 evenly spaced holes with a diameter of 1.27 ± 0.1 cm. The center of each hole is located 2.55 ± 0.03 cm from the center of the disk. 4.5 to 7.0 mm thick Viton sheet, or other suitable autoclavable material, cut into a disk with a diameter of 7.0 ± 0.2 cm containing 6 evenly spaced holes with a diameter of 1.27 ± 0.15 cm (the holes in the Viton are aligned with the holes in the PTFE) and a small hole in the center to house the disk retrieving port. PTFE washer with disk retrieving port. Four nylon screws. PTFE coated 4.0 by 1.4 cm star head magnetic stir bar, set flush against PTFE disk and with holes drilled for assembly with nylon screws. The PTFE ridges on one side of the magnet may be shaved to provide a flush mounting surface. Assemble the pieces conforming to the general details shown in Fig. 2.

NOTE 3—Nominal implies that the manufacturer's tolerance is acceptable.

6.22.4 *Cylindrical Polycarbonate Coupons*, with a diameter of 1.27 ± 0.013 cm and a height of 1.5 to 4.0 mm.

6.23 *Carboys*, two 15 to 20 L autoclavable carboys, to be used for waste and nutrients.

6.23.1 *Carboy Lids*, two carboy lids: one carboy lid with at least 2 barbed fittings to accommodate tubing ID 3.1 mm (one for nutrient line and one for bacterial air vent). One carboy lid with at least 2-1 cm holes bored in the same fashion (one for effluent waste and one for bacterial air vent (filter)).

NOTE 4—Carboy tops can be purchased with fittings.

6.23.2 *Bacterial Air Vent (Filter)*, autoclavable 0.2 μm pore size, to be spliced into tubing on waste carboy, nutrient carboy and reactor top, recommended diameter 37 mm.

7. Reagents and Materials

7.1 *Purity of Water*—All reference to water as diluent or reagent shall mean distilled water or water of equal purity.

7.2 *Culture Media*:

7.2.1 *Bacterial Liquid Growth Broth*, soybean-casein digest medium, or an equivalent general bacterial growth medium. Tryptic Soy Broth is recommended.

⁶ Rotating disk reactor is available commercially from BioSurface Technologies, Corp. www.imt.net/~mitbst, or the user may build the reactor.

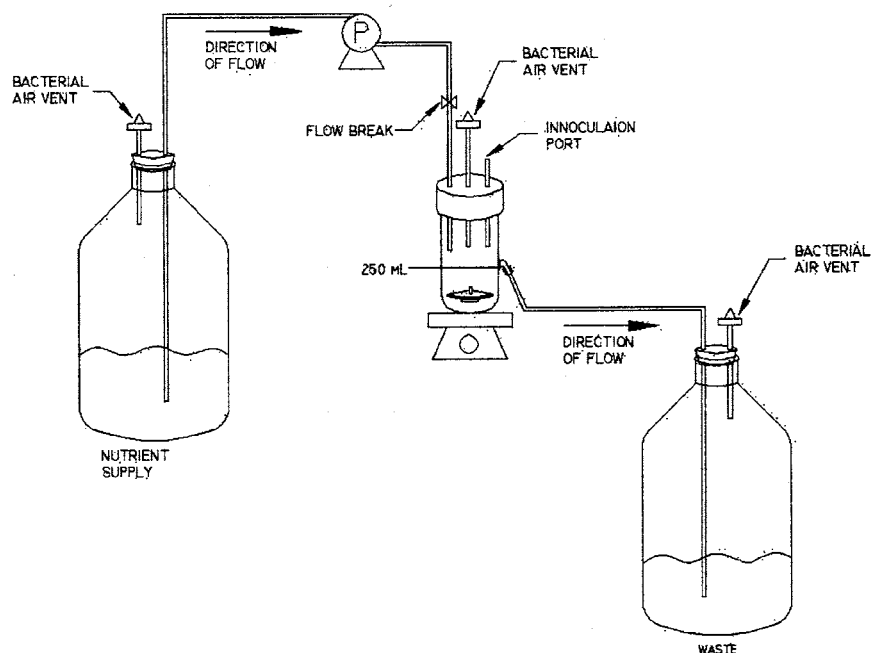


FIG. 1 Rotating Disc Reactor System

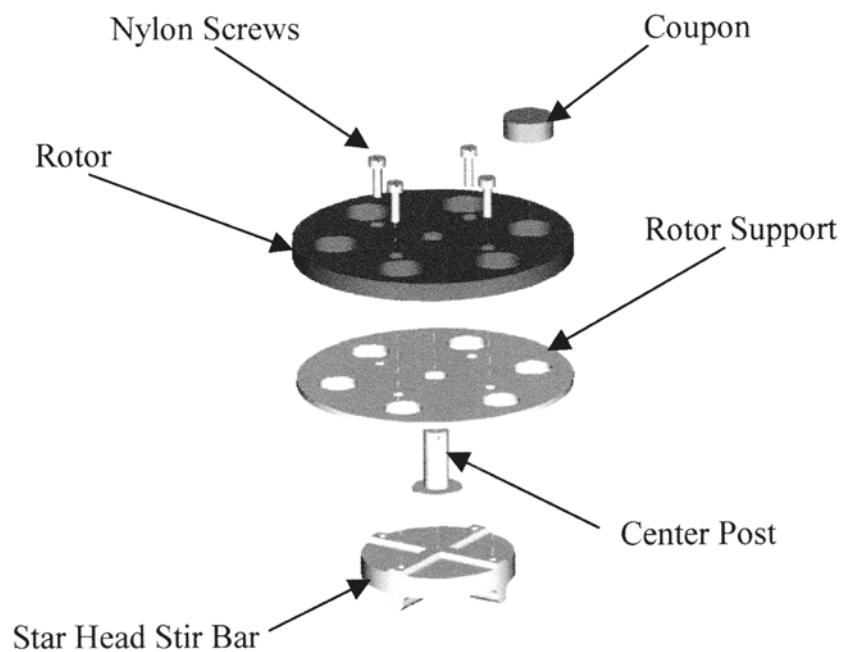


FIG. 2 Rotor Assembly

7.2.2 *Bacterial Plating Medium*, R2A Agar is recommended.

NOTE 5—Media concentration in this protocol differs from the manufacturer's recommendation. Two different concentrations are used in the protocol, 300 mg/L for the inoculum and batch reactor operation and 30

mg/L for the continuous flow reactor operation.

7.3 *Buffered Water*, 0.0425 g/L KH_2PO_4 distilled water, filter sterilized and 0.405 g/L $\text{MgCl} \cdot 6\text{H}_2\text{O}$ distilled water, filter sterilized.

NOTE 6—Prepared according to the document referenced in 2.1.

8. Culture Preparation

8.1 *Pseudomonas aeruginosa* ATCC 700888 is the organism used in this test. An isolated colony is aseptically removed from an R2A plate and placed into 100 mL of sterile bacterial liquid growth broth (300 mg/L) and incubated in an environmental shaker at $35 \pm 2^\circ\text{C}$ for 20 to 24 h. Viable bacterial density should be about 10^8 CFU/mL, which may be checked by serial dilution and plating.

9. Reactor Preparation

9.1 Preparation of Polycarbonate Coupons:

NOTE 7—Coupons can either be used once and discarded or used repeatedly with proper cleaning and sterilization between each use. Check each coupon for scratching, chipping, other damage or accumulated debris before each use by screening under a dissecting microscope at a magnification of at least 20X. Discard those with visible damage to surface topography.

9.1.1 Sonicate coupons for 30 s in a 1:100 dilution of laboratory soap and tap water. The soapy water must completely cover the coupons.

9.1.2 Rinse coupons with reagent water and sonicate for 30 s in reagent water.

9.1.3 Repeat rinsing and sonication with reagent water until, upon visual inspection, no soap is left on the coupons.

9.1.4 Place a coupon into each hole in the rotor, leaving the top of the coupon flush with the Viton surface.

9.1.5 Place the rotor inside of the beaker.

9.2 Preparation of Reactor Stopper Top:

9.2.1 Connect the bacterial air vent by fitting the vent to a small section of appropriately sized tubing, which is then attached to one of the pieces of glass tubing on the reactor stopper top.

9.2.2 The glass flow break is spliced into the nutrient tubing line near the reactor stopper top. See Note 8.

NOTE 8—The other ports on the reactor top may be adapted for anaerobic use, dilution water, or treatment solutions as necessary.

9.3 Sterilization of Reactor System:

9.3.1 Assemble the reactor, with the overflow (waste) line clamped and the reactor top securely fastened to the beaker before sterilization.

9.3.2 Cover the end of the nutrient tubing that connects to the nutrient carboy and the end of the overflow (waste) tubing with aluminum foil. Cover any extra openings on the reactor top with aluminum foil.

9.3.3 Sterilize the reactor system.

10. Procedure

10.1 Batch Phase:

10.1.1 Prepare batch nutrient broth by dissolving bacterial liquid growth medium (300 mg/L) in 250 mL reagent grade water, sterilize.

10.1.2 The top of the sterile reactor is removed and 200 mL of the batch nutrients are aseptically poured into the beaker. The reactor top is secured and the reactor is placed onto a stir plate.

10.1.3 Inoculate the reactor with 1 mL of bacteria from the culture prepared previously (see 8.1). Aseptically inject the

inoculum into the beaker through the inoculation port using the sterile needle and syringe.

10.1.4 The magnetic stir plate is turned on to allow the rotor to spin freely. The rotation speed should equal approximately 200 revolutions per minute (± 100 r/min). The reactor system is allowed to incubate in batch mode at room temperature ($21 \pm 2^\circ\text{C}$) for 24 h.

NOTE 9—The rotational speed of the disk dictates the fluid shear in the reactor. Fluid shear affects biofilm formation.

10.2 Continuous Flow Operation (CSTR Mode):

10.2.1 Aseptically connect the nutrient tubing line to the carboy containing the continuous flow nutrient broth.

10.2.2 Prepare continuous flow nutrient broth by dissolving bacterial liquid growth medium (30 mg/L) in 15 L sterile reagent grade water. Dissolve and sterilize the broth in a smaller volume to prevent caramelization. Aseptically pour the concentrated broth into the carboy of sterile water to make a total of 15 L.

NOTE 10—The concentration of nutrients in the bacterial liquid growth medium will affect the biofilm density that accumulates on the coupons.

10.2.3 A continuous flow of nutrients is pumped into the reactor through a pump set at a flow rate equal to 6.7 ± 0.2 mL/min. Unclamp tubing on the drain spout and place the end into the waste carboy. The drain spout at the 250 mL mark on the beaker allows overflow to occur, maintaining a constant bacterial liquid growth broth concentration of 30 mg/L in the reactor during CSTR (continuously stirred tank reactor) mode.

NOTE 11—Flow rate is calculated by dividing the reactor volume by the residence time. The residence time is 30 min. The reactor volume is approximately 200 mL, (see Note 2). The operator should set an exact flow rate based upon the measured fluid volume in the reactor when the rotor is in place to achieve an exact 30 min residence time. The reactor residence time is a critical parameter that is specific to the bacterial species used during the experiment. In order to select for biofilm growth in the reactor, the residence time must be less than the doubling time for the suspended cells. This will result in the suspended cells washing out of the reactor, leaving only biofilm.

10.2.4 The reactor is operated in CSTR mode for 24 h.

10.3 Sampling the Biofilm:

10.3.1 Prepare sampling materials: vortex, homogenizer, culture tubes, pipettes, empty sterile petri dish, wooden applicator sticks, flame sterilized stainless steel dissecting tools, stainless steel hemostat.

10.3.2 The reactor is opened and the rotor is retrieved using a stainless steel dissecting tool via the retrieving port located in the center post. The rotor is placed into an empty, sterile petri dish.

10.3.3 Holding the rotor by the Viton using a flame-sterilized hemostat, a randomly chosen coupon is removed from the rotor using a stainless steel dissecting tool. Care should be given to not disturb the top coupon surface in any way until it is scraped. If other coupons are to be sampled, replace rotor in reactor to prevent the biofilm from drying.

10.3.4 The coupon is secured by holding it with the flame-sterilized hemostat. The top coupon surface is scraped for approximately 15 s using the flat end of a wooden applicator stick held perpendicular to the coupon surface. The stick is rinsed by stirring it in the dilution buffer in the culture tube.

Repeat the scraping and rinsing process 3 to 4 times, ensuring full coverage of the coupon surface.

10.3.5 The coupon is rinsed by holding the coupon at a 60° angle over the culture tube and pipetting 1 mL of sterile dilution water over the surface of the coupon. The final volume in the culture tube is 10 mL. See **Note 12**.

NOTE 12—Each culture tube originally contains 9 mL of buffered water.

10.4 Analyze the Biofilm Sample:

10.4.1 The scraped biofilm sample is homogenized at $20\,500 \pm 5000$ r/min for 30 s. If more than one biofilm sample is taken, rinse the homogenizer probe by homogenizing a dilution blank for 30 s at the same r/min, homogenize a tube containing 70 % ethanol for 15 s, then remove the probe and let the probe set in the ethanol tube for 1 min. Shake any remaining ethanol off the probe, reattach probe and homogenize a dilution blank for 30 s. Homogenize a second dilution blank. The next sample tube is then ready to be homogenized. Always repeat this cleaning process between samples.

NOTE 13—Homogenizing the sample disaggregates the biofilm clumps to form a homogeneous cell suspension. Improper disaggregation will result in an underestimation of the viable cells present in the sample.

10.4.2 Serially dilute the sample.

NOTE 14—The equation in **10.6** is accurate only if the culture tube the biofilm was scraped into is referred to as the 10^0 dilution.

10.4.3 Plate each dilution in duplicate for colony growth using an accepted plating technique such as spread plating or spiral plating.

10.4.4 Incubate the plates for 17 to 20 h at $35 \pm 2^\circ\text{C}$.

10.5 Cell Enumeration:

10.5.1 Count the appropriate number of colonies according to the plating method used.

10.5.2 Calculate the arithmetic mean of the replicate samples plated.

10.6 The log density for one coupon is calculated as follows:

$$\text{LOG}_{10} (\text{CFU}/\text{cm}^2) = \text{LOG}_{10} [(X/B)(V/A)(D)]$$

where:

X = mean CFU,

B = volume plated,

V = volume scraped into,

A = surface area scraped, and

D = dilution.

10.7 Calculate the overall biofilm accumulation by taking the mean of the log densities calculated in **10.6**.

11. Precision and Bias

11.1 Randomization is used whenever possible to reduce the potential for systematic bias.

11.2 During development, the protocol was repeated 27 times. The internal reproducibility standard deviation for this protocol, when the biofilm density is based on a single coupon, was calculated to equal 0.5. The sources of variability were:

- 40 % attributable to technician-to-technician variability,
- 40 % attributable to experiment-to-experiment within technician variability, and
- 20 % attributable to coupon-to-coupon, within experiment variability.

12. Keywords

12.1 analysis; biofilm; coupon; growth reactor; *Pseudomonas aeruginosa*; reactor; sampling

ASTM International takes no position respecting the validity of any patent rights asserted in connection with any item mentioned in this standard. Users of this standard are expressly advised that determination of the validity of any such patent rights, and the risk of infringement of such rights, are entirely their own responsibility.

This standard is subject to revision at any time by the responsible technical committee and must be reviewed every five years and if not revised, either reapproved or withdrawn. Your comments are invited either for revision of this standard or for additional standards and should be addressed to ASTM International Headquarters. Your comments will receive careful consideration at a meeting of the responsible technical committee, which you may attend. If you feel that your comments have not received a fair hearing you should make your views known to the ASTM Committee on Standards, at the address shown below.

This standard is copyrighted by ASTM International, 100 Barr Harbor Drive, PO Box C700, West Conshohocken, PA 19428-2959, United States. Individual reprints (single or multiple copies) of this standard may be obtained by contacting ASTM at the above address or at 610-832-9585 (phone), 610-832-9555 (fax), or service@astm.org (e-mail); or through the ASTM website (www.astm.org).