



Standard Test Method for In Situ Determination of Network Parameters of Crosslinked Ultra High Molecular Weight Polyethylene (UHMWPE)¹

This standard is issued under the fixed designation F2214; the number immediately following the designation indicates the year of original adoption or, in the case of revision, the year of last revision. A number in parentheses indicates the year of last reapproval. A superscript epsilon (ϵ) indicates an editorial change since the last revision or reapproval.

1. Scope

1.1 This test method describes how the crosslink density, molecular weight between crosslinks, and number of repeat units between crosslinks in ultra-high molecular weight polyethylene (UHMWPE) crosslinked by ionizing radiation or by chemical means can be determined by measuring the swelling ratio of samples immersed in o-xylene. Examples of experimental techniques used to make these measurements are discussed herein.

1.2 The test method reported here measures the change in height of a sample specimen while it is immersed in the solvent. Volumetric swell ratios assume that the sample is crosslinked isotropically, and that the change in dimension will be uniform in all directions. This technique avoids uncertainty induced by solvent evaporation or temperature change.

1.3 The values stated in SI units are to be regarded as standard. No other units of measurement are included in this standard.

1.4 *This standard does not purport to address all of the safety concerns, if any, associated with its use. It is the responsibility of the user of this standard to establish appropriate safety and health practices and determine the applicability of regulatory limitations prior to use.*

2. Referenced Documents

2.1 ASTM Standards:²

D2765 Test Methods for Determination of Gel Content and Swell Ratio of Crosslinked Ethylene Plastics

E691 Practice for Conducting an Interlaboratory Study to Determine the Precision of a Test Method

¹ This test method is under the jurisdiction of ASTM Committee F04 on Medical and Surgical Materials and Devices and is the direct responsibility of Subcommittee F04.15 on Material Test Methods.

Current edition approved May 1, 2008. Published June 2008. Originally approved in 2002. Last previous edition approved in 2002 as F2214 – 02. DOI: 10.1520/F2214-02R08.

² For referenced ASTM standards, visit the ASTM website, www.astm.org, or contact ASTM Customer Service at service@astm.org. For *Annual Book of ASTM Standards* volume information, refer to the standard's Document Summary page on the ASTM website.

3. Terminology

3.1 Definitions of Terms Specific to This Standard:

3.1.1 *crosslink density*, v_d —the theoretical average number of crosslinks per unit volume [mol/dm³].

3.1.2 *molecular weight between crosslinks*, M_c —the theoretical average molecular weight between crosslinks [g/mol].

3.1.3 *swell ratio*, q_s —the ratio of the volume of the sample in an equilibrium swollen state to its volume in the unswollen state.

4. Summary of Test Method

4.1 The height of a cubic specimen is measured, and the specimen is placed in a dry chamber. A selected solvent is chosen according to the Flory network theory and is introduced into the chamber. The chamber is heated to the reference temperature. The sample height is monitored as a function of time until steady state (equilibrium) is achieved. The swell ratio is calculated from the final steady state (equilibrium) height and the initial height.

5. Significance and Use

5.1 This test method is designed to produce data indicative of the degree of crosslinking in ultra high molecular weight polyethylene that has been crosslinked chemically or by ionizing radiation.

5.2 The results are sensitive to the test temperature, solvent, and method used. For the comparison of data between institutions, care must be taken to have the same test conditions and reagents.

5.3 The data can be used for dose uniformity analysis, fundamental research, and quality assurance testing.

6. Apparatus

6.1 The apparatus shall include any device that allows a non-invasive measurement of the change in one dimension of the sample as it swells in the solvent. This measurement could include, but is not limited to:

6.1.1 Mechanical measurements, such as linear variable displacement transducers (LVDTs).

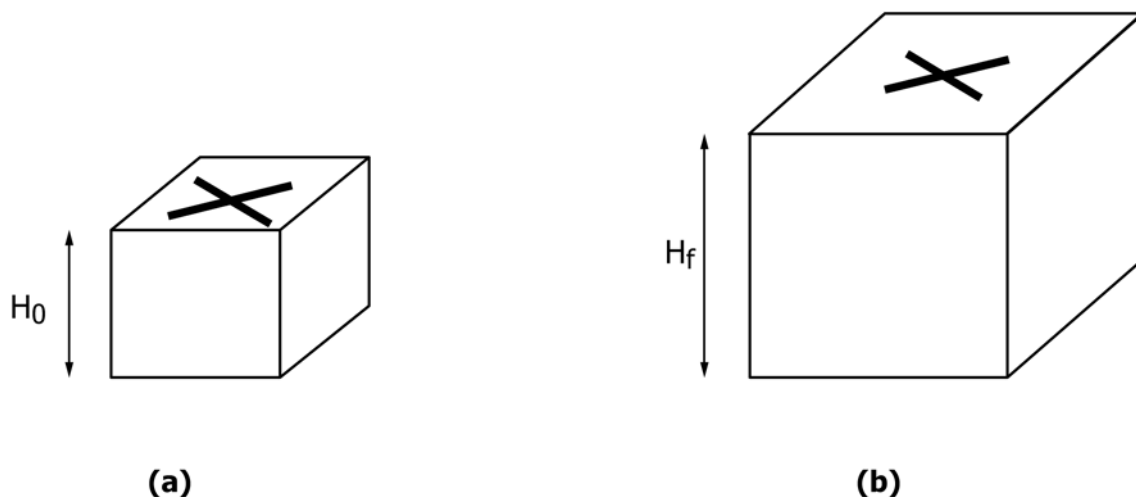


FIG. 1 Marked Measurement Direction Before (a) and After (b) Swelling

6.1.1.1 If a mechanical probe is used, it must be constructed of a material that exhibits little thermal expansion, such as quartz or ceramic.

6.1.2 Optical measurements, such as cameras or laser micrometers.

6.1.2.1 Optical measurements should be insensitive to any refractive index changes in the UHMWPE sample, given the changing temperature of the system.

6.1.3 Inductive measurements, such as proximity sensors. Inductive measurements must be insensitive to temperature or solvent composition.

6.2 The sensitivity of the measurement shall be 1 % of the initial height of the sample, H_0 . An uncertainty analysis has demonstrated that this sensitivity will produce a relative error in crosslink density less than 10 % for samples swollen to a fraction 50 % beyond their initial height. Thicker samples will allow a less sensitive measurement.

6.3 The solvent in the temperature chamber shall be able to reach a temperature of at least 150°C, with an expanded uncertainty of $\pm 1^\circ\text{C}$. Gradients shall not exceed 0.2°C/cm (NB o-xylene boils at 144°C).

6.4 The smallest chamber dimension shall be at least three times the size of the largest initial sample dimension.

6.5 The volume of the chamber shall be at least ten times that of the sample. The chamber should be sufficiently sealed as to prevent gross solvent evaporation during the course of the experiment (typically 2 h).

NOTE 1—The data acquisition software should collect both sample dimension and temperature at a rate of at least 0.1 Hz.

7. Reagents

7.1 *Ortho-Xylene (o-xylene)*, ACS grade, boiling point 144°C.

7.2 *Anti-oxidant*, 2,2'-methylene-bis (4-methyl-6-tertiary butyl phenol).³

8. Safety Precautions

8.1 O-xylene is toxic and flammable, and should be handled only with heat and chemically protective laboratory gloves. The swelling apparatus should ideally be placed inside a vented fume hood, or vented with an elephant trunk should space considerations be an issue. Do not inhale the o-xylene vapors, as dizziness or a headache could result.

8.2 Irganox 1010, the antioxidant, is identified by the manufacturer as an irritant and an inhalation hazard.

9. Test Specimens

9.1 At least three specimens with a minimum sample height of 500 μm should be machined. The top and bottom surfaces should be parallel and smooth. The width and length (or diameter, in the case of cylindrical samples) should be less than one third the size of the sample chamber (see 6.4). The height-to-width aspect ratio should be at least 1:2 to minimize buckling, with 1:1 preferred. The machining should be performed so as to minimize thermal degradation of the samples.

9.2 *Orientation of Samples*—Given that the swelling behavior can depend on molecular alignment induced by processing conditions, the test specimens should be machined so that the relevant processing direction can be easily identified. The samples can then be oriented in the swelling apparatus relative to the molding direction (that is, perpendicular to the extrusion of compression molding direction). The specimens can be marked as shown in Fig. 1 to aid in sample alignment.

10. Procedure

10.1 Add approximately 0.5 to 1 % (mass fraction) of the antioxidant to the o-xylene to make a stock solution.

³ Trade name: Irganox 1010 has been found satisfactory for this purpose. Available from Ciba-Geigy, 540 White Plains Rd., P.O. Box 2005, Tarrytown, NY 10591-9005.

10.2 The initial sample height should be measured with a resolution of 1 % of the sample height using a micrometer. This value should be recorded. The measurement direction on the sample can be indicated with a permanent marker. An example is shown in Fig. 1.

10.3 The sample should be pre-wet with o-xylene, then quickly placed in the dry chamber with the sample correctly oriented as marked in 10.2.

10.4 The initial sample dimension, as determined with the measurement system of the instrument described in 6.1, should be recorded.

10.5 Start recording the sample dimension at a minimum rate of 1 point every 10 s.

10.6 Introduce the o-xylene stock solution into the chamber at a slow rate to prevent disturbing the sample.

10.7 Raise the temperature of the solvent in the chamber to $130 \pm 1^\circ\text{C}$.

10.8 Continue to monitor the temperature and sample dimension until equilibrium is achieved (within $\pm 10 \mu\text{m}$) over a period of 15 min.

10.9 Decrease the temperature to below 50°C . Discard the o-xylene in an environmentally responsible manner, and clean the sample cell thoroughly.

10.10 Examine the sample after the test is complete. If it has shown signs of cracking, or is yellowed, thermal degradation is likely to have occurred. This data will be suspect and should be discarded.

11. Calculation of Swell Ratio

11.1 The swell ratio, q_s , is computed as indicated from the height measurement:

$$q_s = (V_f/V_0) = (H_f/H_0)^3 \quad (1)$$

where:

V_f = final volume,
 V_0 = initial volume,
 H_f = final height, and
 H_0 = initial height.

NOTE 2—This calculation assumes that the sample is isotropic.

12. Calculation of Crosslink Density and Molecular Weight Between Crosslinks

12.1 Given the steady state swell ratio, q_s , of a polymer immersed in a specific solvent at a particular temperature, the crosslink density, molecular weight between crosslinks, and number of crosslinks/chain can be computed if one knows the Flory interaction parameter, χ_1 , for the polymer-solvent system.

12.2 From Flory's network theory, which explains the swell ratio of a polymer-solvent system as a competition between elastic forces and forces derived from the free energy of mixing, the following expression is derived for the crosslink

TABLE 1 Summary of Mean ($\pm S_R$, Absolute Interlaboratory Uncertainty) Swell Ratio (q), Crosslink Density (v_d), and Molecular Weight Between Crosslinks (M_c) for the Four Samples

Dose, kGy	Swell Ratio, q	Crosslink Density, v_d [mol/dm ³]	Molecular Weight Between Crosslinks, M_c [g/mol]
54.2	3.37 ± 0.26	0.133 ± 0.017	$7,650 \pm 1,010$
71.5	3.12 ± 0.24	0.151 ± 0.021	$6,720 \pm 920$
89.2	3.12 ± 0.24	0.152 ± 0.020	$6,700 \pm 890$
110.1	2.97 ± 0.31	0.170 ± 0.032	$6,150 \pm 1,190$

density, v_d , as a function of the steady state swelling ratio, the Flory interaction parameter, and ϕ_1 , the molar volume of the solvent.^{4,5}

$$v_d = -\frac{1n(1 - q_s^{-1}) + q_s^{-1} + \chi_1 q_s^{-2}}{\phi_1(q_s^{-1/3} - q_s^{-1/2})} \quad (2)$$

12.2.1 The expression in Eq 2 assumes a three-dimensional network composed primarily of "H-bonding," or the formation of crosslinks along the main chain rather than at the chain ends. Additionally, network entanglements may partially contribute to the refractive forces. Thus the calculated crosslink density will account for these contributions as well. The expression in Eq 2 has been shown to be valid for swelling ratios up to $q = 10$, or $M_c > 10\,000$ g/mol.

12.3 The expression in Eq 2 can be reduced to calculate the molecular weight between crosslinks, M_c , where $\bar{v}$ is the specific volume of the polymer.

$$M_c = (\bar{v}v_d)^{-1} \quad (3)$$

12.4 For polyethylene in o-xylene at 130°C , the following approximate parameters shall be used:

Parameter	Value
χ_1 (o-xylene-PE, 130°C)	$0.33 + 0.55/q_s$
ϕ_1 (o-xylene) [cm ³ /mol]	136
$\bar{v}^{-1}$ [g/dm ³]	920

13. Report

13.1 Report the following information:

- 13.1.1 Complete identification of the sample,
- 13.1.2 Solvent and temperature used,
- 13.1.3 Initial heights of the samples,
- 13.1.4 Final heights of the samples,
- 13.1.5 Calculated swell ratio, crosslink density, and molecular weight between crosslinks, and
- 13.1.6 The orientation of the specimens relative to the principle processing direction (that is, ram extrusion direction, or compression molding direction).

14. Precision and Bias

14.1 Tables 1 and 2 are based on a round robin study

⁴ Flory, P. J., *Principles of Polymer Chemistry*, Ithaca and London, Cornell University Press, 1953.

⁵ Flory, P. J., and Rehner, J., "Statistical mechanics of cross-linked polymer networks. II. Swelling," *J. Chem. Phys.*, Vol 11, No. 11, 1943, pp. 521–526.

TABLE 2 Summary of Standard Interinstitutional (s_R) and Intrainstitutional (s_r) Relative Uncertainty for Swell Ratio (q), Crosslink Density (v), and Molecular Weight Between Crosslinks (M_c)

Dose, kGy	Uncertainty, q		Uncertainty, v_d		Uncertainty, M_c	
	s_r (%)	s_R (%)	s_r (%)	s_R (%)	s_r (%)	s_R (%)
54.2	3.0	7.8	4.7	12.8	5.1	13.2
71.5	2.3	7.8	4.3	13.6	4.0	13.7
89.2	1.9	7.5	3.7	13.5	3.3	13.3
110.1	5.4	10.5	9.4	18.7	10.1	19.3

conducted in 2001 involving four sets of ultra high molecular weight polyethylene test samples tested by six laboratories.⁶ For all sets of samples, all the specimens were prepared at the same time by the same laboratory. Each laboratory tested three specimens from each set of samples.

14.1.1 Samples from NIST Reference® 8456, Ultra High Molecular Weight Polyethylene were machined into 5 mm cubes, packaged in nitrogen, and irradiated four gamma different irradiation doses (54.2, 71.5, 89.2, and 110.1 kGy).⁷

⁶ Spiegelberg, S., Kurtz, S., Muratoglu, O., Greer, K., Costa, L., Wallace, S., and Cooper, C., *Interlaboratory Reproducibility of Swell Ratio Measurements for Crosslinked Polyethylene*, 48th Annual Meeting of the Orthopedic Research Society, Dallas, TX, 2002.

⁷ Cubes of 5 mm per side are available from NIST as Reference® 8457 Orthopaedic Grade Polyethylene Cubes. Available from National Institute of Standards and Technology (NIST), 100 Bureau Dr., Stop 3460, Gaithersburg, MD 20899-3460.

14.2 The data indicates that for the range of irradiation doses currently used to produce highly crosslinked UHMWPE for orthopedic applications, the swell ratio measurement is associated with interlaboratory standard uncertainty of 8 to 11 %.

14.3 *Concept of r and R* —If S_r and S_R , the absolute intralaboratory and interlaboratory uncertainties, have been calculated from a large enough body of data, and for test results that are averages from testing 3 specimens:

14.3.1 *Repeatability, r* —Comparing two test results for the same material, obtained by the same operator using the same equipment on the same day, the two test results should be judged not equivalent if they differ by more than the r value for that material.

14.3.2 *Reproducibility, R* —Comparing two test results for the same material, obtained by different operators using different equipment on different days, the two test results should be judged not equivalent if they differ by more than the R value for that material.

14.3.3 Any judgment in accordance with 14.3.1 and 14.3.2 would have an approximate 95 % (0.95) probability of being correct.

14.4 There is no recognized standard by which to estimate bias of this test method.

15. Keywords

15.1 crosslink density; molecular weight; network; swelling; UHMWPE

APPENDIXES

(Nonmandatory Information)

X1. RATIONALE

X1.1 This test method is intended to describe the procedure to be followed in order to measure the swelling behavior of crosslinked ultra high molecular weight polyethylene.

X1.2 There are contributions to the measured swelling behavior introduced when the polyethylene melts, as determined by the change in density as the sample goes from a solid to a melt. The relative or fractional change in height from this contribution has been calculated to be 2.5 %.⁸

X1.3 There are contributions from thermal expansion of the sample, which can be calculated from the thermal coefficient of expansion of UHMWPE. The relative change over 130°C is less than 0.05 %.

X1.4 The weight of the sample probe, for contact measurements, can cause erroneously low measured swelling

ratios if the probe constrains the sample. This contribution will depend on the modulus of the sample, which in turn depends on the degree of swelling. As an example, a sample with a modulus of 100 MPa, when subjected to a pressure of 1000 Pa, will undergo a decrease in the height of 1 %.

X1.5 Polymer molecules can align in flow fields introduced during normal processing conditions, such as compression molding or extrusion. The molecular alignment can be frozen into the consolidated polymer, yielding a lower entropic state than a fully randomly oriented polymer system would exhibit. If a test specimen is machined from a block of anisotropically oriented material and subjected to a swelling test, the polymer chains in the aligned direction will deform more than the less-aligned direction, yielding a different linear swelling behavior in these orthogonal directions. Consequently, to improve repeatability/reproducibility in tests, the test orientation relative to the principle processing directions must be the same for all test specimens. Users can measure the swell ratio in either direction, and report both results.

⁸ Muratoglu, O. K., Cook, J. L., et al., *A Novel Technique to Measure the Crosslink Density of Irradiated UHMWPE*, 24th Annual Meeting of the Society for Biomaterials, San Diego, CA, 1998.

X2. INDEX TO SELECTED STATISTICAL TERMS (per Practice E691)

s_r —Repeatability standard deviation
 s_R —Reproducibility standard deviation
 S_r —Absolute repeatability uncertainty
 S_R —Absolute reproducibility uncertainty

ASTM International takes no position respecting the validity of any patent rights asserted in connection with any item mentioned in this standard. Users of this standard are expressly advised that determination of the validity of any such patent rights, and the risk of infringement of such rights, are entirely their own responsibility.

This standard is subject to revision at any time by the responsible technical committee and must be reviewed every five years and if not revised, either reapproved or withdrawn. Your comments are invited either for revision of this standard or for additional standards and should be addressed to ASTM International Headquarters. Your comments will receive careful consideration at a meeting of the responsible technical committee, which you may attend. If you feel that your comments have not received a fair hearing you should make your views known to the ASTM Committee on Standards, at the address shown below.

This standard is copyrighted by ASTM International, 100 Barr Harbor Drive, PO Box C700, West Conshohocken, PA 19428-2959, United States. Individual reprints (single or multiple copies) of this standard may be obtained by contacting ASTM at the above address or at 610-832-9585 (phone), 610-832-9555 (fax), or service@astm.org (e-mail); or through the ASTM website (www.astm.org). Permission rights to photocopy the standard may also be secured from the Copyright Clearance Center, 222 Rosewood Drive, Danvers, MA 01923, Tel: (978) 646-2600; <http://www.copyright.com/>